Panasonic

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Panasonic

Operating Instructions (Household) Sonic Vibration Toothbrush Model No. EW-DL34

Contents

Safety precautions	4	Changing the brush head (Sold separately)	15
Parts identification	8	FAQ	15
Before cleaning your teeth	8	Troubleshooting	16
Features	8	Battery life	18
Charging the toothbrush	9	Removing the built-in rechargeable battery	18
Basic usage	10	Specifications	18
Tips on usage	12		
After use	13		

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

DANGER

Denotes a potential hazard that will result in serious injury or death.

WARNING

Denotes a potential hazard that could result in serious injury or death.

CAUTION

Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe the type of instructions to be observed.

This symbol is used to alert users to a specific operating procedure that must not be performed.

This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

WARNING

► Power supply

Do not connect or disconnect the power plug to a household outlet with a wet hand.

- Doing so may cause electric shock or injury.

Do not immerse the charger in water or wash it with water.

Do not place the charger over or near water filled sink or bathtub.

Do not use when the cord or the power plug is damaged or when the fitting into the household outlet is loose.

Do not damage the cord or the power plug.

• **Do not deface, modify, forcefully bend, pull, or bundle the cord or the power plug. Also, do not place anything heavy on, twist or pinch the cord.**

- Doing so may cause electric shock or fire due to a short circuit.

Do not use anything other than the supplied charger for any purpose.

Also, do not use any other product with the supplied charger. (See page 8.)

- Doing so may cause burn or fire due to a short circuit.

Do not use in a way exceeding the rating of the household outlet or the wiring.

- Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.

WARNING

Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the charger.

Fully insert the power plug.

- Failure to do so may cause fire or electric shock.

Always unplug the power plug from a household outlet when cleaning the charger.

- Failure to do so may cause electric shock or injury.

Regularly clean the power plug to prevent dust from accumulating.

- Failure to do so may cause fire due to insulation failure caused by humidity.

Disconnect the power plug and wipe with a dry cloth.

► In case of an abnormality or malfunction

Immediately stop using and remove the power plug if there is an abnormality or failure.

- Failure to do so may cause fire, electric shock, or injury.

<Abnormality or failure cases>

• **The main unit or charger is deformed or abnormally hot.**

• **The main unit or charger smells of burning.**

• **There is abnormal sound during use or charging of the main unit or charger.**

- Immediately request inspection or repair at an authorized service center.

► This product

This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, or leave in a high temperature environment.

- Doing so may cause overheating, ignition, or explosion.

Do not modify or repair.

- Doing so may cause fire, electric shock, or injury. Contact an authorized service center for repair (battery change, etc.).

Never disassemble except when disposing of the product.

- Doing so may cause fire, electric shock, or injury.

► Preventing accidents

Do not store within the reach of children or infants. Do not let them use it.

- Doing so may cause an accident or injury due to accidental ingestion of removable parts.

► Take care of the following precautions as well

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

- Failure to do so may cause an accident or injury.

WARNING

-  **The supply cord cannot be replaced. If the cord is damaged, the charger should be scrapped.**
- Failure to do so may cause an accident or injury.

CAUTION

▶ Preventing damage to teeth and gums

- Do not use the brush with toothpaste that contains a high abrasive content, such as toothpaste for removing tar or for whitening. Do not apply the brush head strongly against teeth or gums, or use it continuously against same spot.**
-  **Do not apply the plastic part of the main unit or brush head to teeth.**
Do not let anyone who cannot operate the brush by themselves use it.
People with low sensitivity in their mouth should not use the brush.
- Doing so may cause injury or damage to teeth and gums.

- People who may have periodontal disease, have teeth being treated, or who are concerned about symptoms in their mouth should consult a dentist before use.**
- Failure to do so may cause injury or damage to teeth and gums.

▶ This product

-  **Never use for a purpose other than cleaning your teeth.**
- Doing so may cause an accident or injury.
-  **Do not share the brush head with your family or other people.**
- Doing so may result in infection or inflammation.
- ▶ **Note the following precautions**
 -  **Do not drop or hit against other objects.**
- Doing so may cause injury.
 - Do not wrap the cord around the charger when storing.**
 - Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.
 - Disconnect the power plug from the household outlet when not charging.**
 - Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.
 - Disconnect the power plug by holding onto it instead of the cord.**
 - Failure to do so may cause electric shock or injury.

DANGER

- The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.
 - Do not charge the battery after it has been removed from the product.
 - Do not throw into fire or apply heat.
 - Do not hit, disassemble, modify, or puncture with a nail.
 - Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
 - Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
 - Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
 - **Never peel off the tube.**
- Doing so may cause overheating, ignition, or explosion.

WARNING

- After removing the rechargeable battery, do not keep it within reach of children and infants.**
- The battery will harm the body if accidentally swallowed.
If this happens, consult a physician immediately.

- If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.**
- Failure to do so may cause eye injury.
Consult a physician immediately.

CAUTION


- If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.**
- Failure to do so may cause inflammation.

Note

- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type references of power supply unit is marked near the symbol.

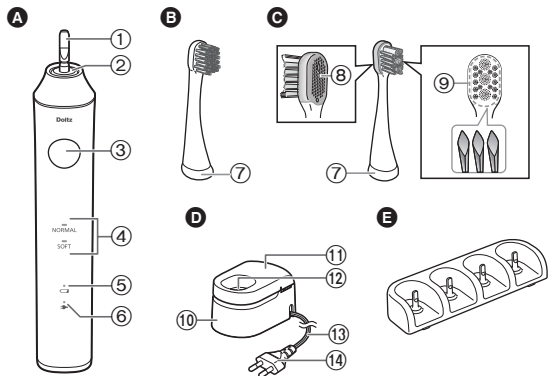Parts identification

A Main unit

- ① Toothbrush shaft
- ② Waterproof rubber seal
- ③ Power switch
(power off/on/mode selection)
- ④ Mode indicator
- ⑤ Battery capacity lamp
()
- ⑥ Charge status lamp
()

B Compact extra fine brush

- ⑦ Identification ring

C Compact stain care brush

- ⑧ Tongue cleaner
- ⑨ High-density bristles

D Charger (RE8-56) (The shape of the power plug differs depending on the area.)

- ⑩ Charger base
- ⑪ Main unit guide
- ⑫ Charging section
- ⑬ Cord
- ⑭ Power plug

Accessory

E Brush stand

Before cleaning your teeth

► This product

The Sonic Vibration Toothbrush is an electric toothbrush with sonic vibration* of approximately 31 000 horizontal brush strokes/minute.

* Sonic vibration in acoustic wave field

► Bleeding of the gums

When you first begin using the toothbrush, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the toothbrush for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the toothbrush and consult your dentist.

Features

► About soft start function

To get accustomed to the high speed vibration, the toothbrush will gradually reach the selected speed and vibration approximately 2 seconds after it starts operating.

► About the 30-second brushing notification timer function

The toothbrush will pause once at 30-second intervals to indicate the brushing time. It pauses twice at 2-minute intervals then starts again.

► About mode storage function

The mode used last is stored.

- The toothbrush operates in the mode last used when the power switch is turned on.


Charging the toothbrush

Preparation

Place the charger on a horizontal and stable place.

1 With the main unit facing forwards, set straight into the charger.


- Charging cannot start when the main unit is tilted or the bottom of the main unit is not in contact with the charging section.

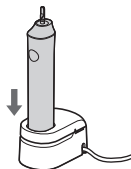► Charging timing guidelines


- It is time to charge when the battery capacity lamp (☐) starts to blink while using.
- If you set the switch to "OFF" with low remaining charge, the battery capacity lamp flashes and the alarm beeps.
- The battery capacity lamp turns off after about 5 seconds.

Charging reminder

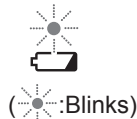2 Insert the power plug into the household outlet.

- Charging is completed after approx. 17 hours.
(The charging time may differ with the charging capacity.)

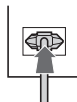► Display while charging

While charging	After charging is completed	Abnormal charging
● ☐	○ ☐	☐ ☐
		The charge status lamp (☐) will blink twice every second.

(●:Glow/☐:Blinks/○:Turns off)

3 Disconnect the power plug from the household outlet.

(for safety and to reduce energy usage)

- Charging the battery for more than 17 hours will not damage the main unit.

► Battery runtime on a full charge

Mode	Battery runtime
NORMAL SOFT	Approx. 90 minutes

(Use time may vary depending on ambient temperature and use conditions.)

Notes

- If the charge status lamp (☐) blinks quickly (twice per second), it may be due to a charging error or because the main unit is not properly mounted on the charger. Confirm the recommended temperature and mount the main unit on the charger properly.

- Recommended ambient temperature for charging is 5 °C – 35 °C. The battery may not charge properly or not at all under extreme low or high temperatures.
- The charging time may vary when using for the first time, or when it has not been in use for more than 6 months. Also, the charge status lamp (➡) may not glow for a few minutes after charging starts, but it will glow if charging continues.


Basic usage

Preparation


Select the appropriate brush for the desired use. (See page 12.)

1 Attach the brush head to the main unit.

- Align the front of the main unit and the front of the brush head.
1. Ensure that the direction is correct.
 - Do not attach the brush head by twisting it into place. Doing so may damage the brush head and main unit.
 - The toothbrush shaft will move when a force is applied. This is required for operation, and is not a malfunction.
 2. Apply a small amount of toothpaste.

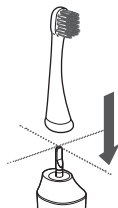There is a small gap between the brush head and the main unit (a few mm), but this is not a malfunction.

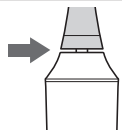2 To avoid spattering of toothpaste, put the brush head in your mouth and press the power switch to turn on the main unit.

► Select the desired mode

- The device starts in the previous mode used when the power is turned on.
- Press the switch within 5 seconds to change modes. (The blue mode indicator illuminates to indicate the selected mode.)

To switch the mode, press the power switch within 5 seconds.

NORMAL ←
 ▼
 SOFT OFF


3 Lightly press and move the toothbrush against your teeth.

► Example brushing sequence

Divide your teeth into 4 parts (upper and lower, left and right), and brush one part in order such as “① front of the teeth”, “② biting plane of the teeth”, and “③ back of the teeth” spending 30 seconds.

Repeat this same procedure for the other 3 areas.

- By brushing the above 4 parts for 30 seconds each (about 2 minutes), you can evenly brush your teeth without missing any spots in your mouth.
- See page 8 about “30-second brushing notification timer function.”

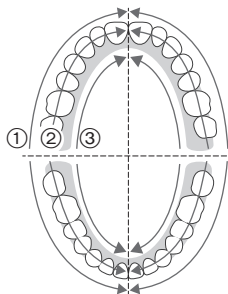4 Turn off the main unit before removing it from your mouth to prevent saliva from spattering.

- Press the switch after at least 5 seconds of operation in one mode to turn the toothbrush off. If you press the switch before 5 seconds of operation have elapsed, the mode changes. (See page 10.)

Notes

- Do not put the main unit in contact with hair products, hand creams, etc. Doing so may damage the exterior.
- The appliance may not operate in an ambient temperature approximately 5 °C or lower.
- Use SOFT mode if it feels too strong in NORMAL mode.
- When first using this appliance, some people may experience a tickling or tingling sensation caused by the sonic vibrations. The sensation should subside after you use the toothbrush several more times.
- If you push the brush head too hard, the sound may change, but this is not an abnormality.

Tips on usage

Change the orientation and manner in which you brush your teeth in accordance with your objective and the area you are brushing.

Periodontal pocket care

Recommended mode: NORMAL

Compact extra fine brush

When you want to brush mainly the periodontal pockets or between the teeth and gums.

- You can also brush the teeth front surface and the biting surface.

Stain care

Recommended mode: NORMAL

Compact stain care brush

When you want stain care during normal brushing.
(This does not mean that all stains can be removed.)

How to apply/move the toothbrush

Apply the bristles of the brush lightly against the teeth and move the toothbrush slowly.

Front surface and the biting surface of teeth

- Apply the brush head at a 90° angle.

Periodontal pocket and the boundary between the teeth and gums

- Use the compact extra fine brush.
- Place the brush head at a 45° angle.
- Apply the bristles lightly into the periodontal pocket and move the toothbrush slowly.

For the back of the front teeth or misaligned teeth

- Place the brush vertically.
- Move toothbrush to match the angle of the teeth.

Tongue care

Recommended mode: SOFT

Tongue cleaner

Compact stain care brush (Back)

When you want to clean the tongue.

How to apply/move the toothbrush

Apply the tongue cleaner so that it barely touches the tongue, and move forward with slight force from the area that becomes dirty.

Tongue

- Repeat several times.
- Poke your tongue out forward and apply the brush while you are breathing out to prevent emetic response.

Note

- You may injure your tongue if the tongue cleaner is pressed against the tongue too strong.

After use

After using the toothbrush, clean the toothbrush to maintain cleanliness of the device.

- Use neutral kitchen detergent when cleaning the toothbrush with detergent.
- Do not use thinner, benzene, alcohol, etc. to clean the appliance.
Doing so may cause failure, cracking or discoloration of the parts.

Remove the brush head.

Remove the brush head straight from the toothbrush shaft.

Cleaning the brush head

Wash the brush head under running water.

Wash the soiling on the inside and bottom areas of the brush head thoroughly.


- Do not forcefully pull the bristles of the brush or bend them out excessively.
- Do not wash the tongue cleaner excessively with pointed objects or your fingernails.
- Do not wash with water hotter than 80 °C. Doing so may damage the brush head.

Cleaning the main unit

Wash the main unit under running water.

- Do not damage or rub the waterproof rubber seal (A) with hard materials.
Do not leave the main unit under water. Water may enter inside.
- Do not wash with water hotter than 40 °C. Doing so may damage the main unit.

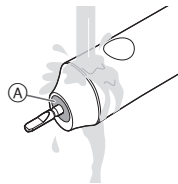Cleaning the charger

1. Disconnect the power plug from the household outlet.
2. Slide and remove the main unit guide from the charger base.


3. Wash the main unit guide with water.

Wipe off water after washing.


- Do not rinse with water except the main unit guide.

4. Wipe off any dirt on the charger base using a soft cloth.

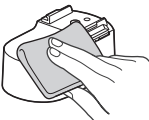► Attaching the main unit guide

Align the main unit guide with the designated lines on the charger base and then slide the main unit guide until it clicks into place.

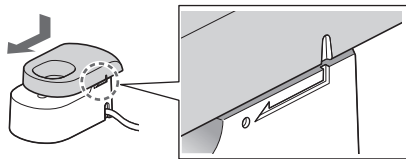Wipe the prongs of the power plug with a dry cloth about once every 6 months.

Cleaning the brush stand

Wash and wipe off water.


► After cleaning

Wipe off water from the brush head and main unit for storage.

- Do not charge while it is wet. Doing so may cause soiling.

► If the identification ring comes off

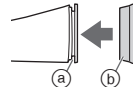Attach by aligning the indentation (a) in the brush head with the raised part (b) of the inside of the identification ring.

- Do not forcefully pull the identification ring. Doing so may cause damage.
- Using the toothbrush without the identification ring has no effect on its performance.

Changing the brush head (Sold separately)

- Brushes are consumables.
- The brush should be replaced with a new one every 3 months or so due to hygiene reasons, even if its shape is unchanged.
Even within 3 months, please exchange the brush head if the bristles are frayed.
- Colored identification rings are included with the separately sold brush.
Optional brush heads with identification ring to allow several persons to share the main unit.

Part number	Replacement Brush for EW-DL34	
WEW0914	Compact extra fine brush	
WEW0938	Compact stain care brush	

FAQ

	Question	Answer
While using	Why does the toothbrush temporarily stop every 30 seconds? ▶	The toothbrush will notify you of the elapsed brushing time every 30 seconds. (See page 8 about “30-second brushing notification timer function.”)
Battery	Will the battery deteriorate if not used for an extended period of time? ▶	The battery will deteriorate if not used for over 6 months, so make sure to perform a full charge at least once every 6 months.
	Can I charge the appliance after every use? ▶	The battery used is a lithium-ion, so charging after every use will have no effect on battery life.

	Problem	Possible cause	Action
While charging	The charge status lamp (➡) does not glow when charging starts.	▶ You have just purchased the toothbrush or have not used it for a while.	▶ The charge status lamp (➡) may not glow for several minutes after starting to charge it, but it will glow if you continue to charge it. If the charge status lamp (➡) does not glow after a while, consult with an authorized service center for repair.
		▶ The main unit is not properly mounted on the charger.	▶ Insert the main unit so the bottom of it contacts the charging section of the charger, and confirm that the charge status lamp (➡) glows. (See page 9.)
	The charge status lamp (➡) is blinking rapidly. (Twice every second)	▶ The temperature is too low or high while charging.	▶ Charge the toothbrush within the recommended charging temperature of 5 °C to 35 °C. ▶ If the charge status lamp (➡) blinks even when charging within the recommended temperature range, please contact an authorized service center for repair.

Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge.

(The operating life of the battery may differ depending on usage or storage conditions.)

Removing the built-in rechargeable battery

Remove the built-in rechargeable battery before disposing the toothbrush.


The battery is to be disposed of safely.

Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing the toothbrush, and must not be used to repair it. If you dismantle the toothbrush yourself, it will no longer be waterproof, which may cause it to malfunction.

- Press the power switch to turn on the power and then keep the power on until the battery is completely discharged.

- Follow the procedure below, using a screwdriver to disassemble the appliance. Perform steps ① to ④ and lift the battery, and then remove it.
- Take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.

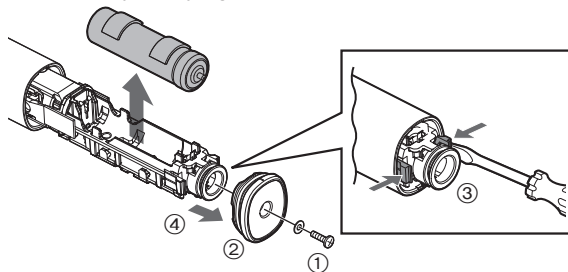For environmental protection and recycling of materials

This toothbrush contains a Li-ion battery.

Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

Power supply	See the name plate on the product (Automatic voltage conversion)
Charging time	Approx. 17 hours (Room temperature: 5 °C – 35 °C)

This product is intended for household use only.

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© Panasonic Corporation 2019
Panasonic Corporation